

# QDrone2

## Power User Manual

*v 1.0 – 9th Jan 2023*

For more information on the solutions Quanser Inc. offers,  
please visit the web site at: <http://www.quanser.com>

Quanser Inc. info@quanser.com  
119 Spy Court Phone : 19059403575  
Markham, Ontario Fax : 19059403576  
L3R 5H6, Canada printed in Markham, Ontario.

This document and the software described in it are provided subject to a license agreement. Neither the software nor this document may be used or copied except as specified under the terms of that license agreement. Quanser Inc. grants the following rights: a) The right to reproduce the work, to incorporate the work into one or more collections, and to reproduce the work as incorporated in the collections, b) to create and reproduce adaptations provided reasonable steps are taken to clearly identify the changes that were made to the original work, c) to distribute and publicly perform the work including as incorporated in collections, and d) to distribute and publicly perform adaptations. The above rights may be exercised in all media and formats whether now known or hereafter devised. These rights are granted subject to and limited by the following restrictions: a) You may not exercise any of the rights granted to You in above in any manner that is primarily intended for or directed toward commercial advantage or private monetary compensation, and b) You must keep intact all copyright notices for the Work and provide the name Quanser Inc. for attribution. These restrictions may not be waved without express prior written permission of Quanser Inc.

**FCC Notice** This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Industry Canada Notice** This Class A digital apparatus complies with Canadian ICES-003. Cet appareil numérique de la classe A est conforme à la norme NMB-003 du Canada.

#### Waste Electrical and Electronic Equipment (WEEE)



This symbol indicates that waste products must be disposed of separately from municipal household waste, according to Directive 2002/96/EC of the European Parliament and the Council on waste electrical and electronic equipment (WEEE). All products at the end of their life cycle must be sent to a WEEE collection and recycling center. Proper WEEE disposal reduces the environmental impact and the risk to human health due to potentially hazardous substances used in such equipment. Your cooperation in proper WEEE disposal will contribute to the effective usage of natural resources.

This product meets the essential requirements of applicable European Directives as follows:

**CE Compliance** 

- 2006/95/EC; Low-Voltage Directive (safety)
- 2004/108/EC; Electromagnetic Compatibility Directive (EMC)

**Warning:** This is a Class A product. In a domestic environment this product may cause radio interference, in which case the user may be required to take adequate measures.



**This equipment is designed to be used for educational and research purposes and is not intended for use by the public.** The user is responsible to ensure that the equipment will be used by technically qualified personnel only. While the end-effector board provides connections for external user devices, users are responsible for certifying any modifications or additions they make to the default configuration.

## Table of Contents

|                                  |   |
|----------------------------------|---|
| A. Charging the battery          | 3 |
| B. Connecting the battery        | 4 |
| C. Using Static Power            | 4 |
| D. Turning ON the QCar           | 5 |
| E. Turning OFF the QCar          | 6 |
| F. Low-battery and auto-shutdown | 6 |

## A. Charging the battery

1. Power the battery charger (Figure 1a) with the supplied power cable (Figure 1b).
2. Connect the balancer cable on the battery to the charger's middle port (Figure 1c).
3. Charging will start automatically. Charging is complete when the fourth LED is illuminated. A full charge takes approximately 40 minutes.
4. After charging is complete, unplug AC power first. Then disconnect the battery from the charger.



a. charger/balancer



b. Power cable



c. connections to charger

Figure 1. Wiring and using the RC battery charger/balancer

## B. Connecting the battery

Once the LiPo battery is charged, place it in the QCar by following these steps.

1. As shown in Figure 2, press down on the locking arm (1) and pull outward on the latch (2).
2. Place the LiPo battery (3) in the compartment with the battery cables towards the front of the QCar. Do not connect the battery cables yet.
3. Push the latch (2) back in over the battery and its cables, and into the lock until the locking arm (1) clicks.
4. Connect the battery's female XT-60 connector into the QCar's male XT-60 connector (4) in Figure 2.

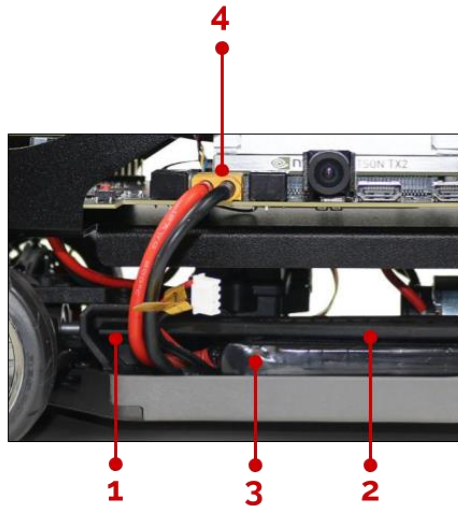


Figure 2. Battery compartment and connectors

## C. Using Static Power

For applications where the QCar platform will remain static, it may be powered using a power supply instead of the provided batteries. To do this, connect the male 6-pin connector of the supplied power supply (Figure 3a) to the female 6-pin connector on the supplied adapter in Figure 3b. Connect the adapter's female XT-60 connector to the QCar's male XT-60 battery connector, and continue operation as normal.



a. Power supply



b. XT-60 adaptor PCB

Figure 3. Static power supply and adaptor


## D. Turning ON the QCar

After the battery is connected, turn on the QCar platform by firmly but gently pressing on the red power button shown in Figure 4.



Figure 4. Red power button and green power LEDs on the QCar PCB

The power LEDs (Figure 4, labelled VBAT, 5V and 3.3V on the PCB) should turn ON. The headlamps, brake lamps, reverse lights and indicators should flash once and the speaker will emit a click sound. The LCD will display 'QUANSER INNOVATE EDUCATE' as shown in Figure 5a. After a few seconds, the LCD will change to 'qcar-xxxxx starting up ...' as shown in Figure 5b. Once the car boots up completely, it will switch to a default mode which shows the current battery level, Ethernet/WiFi status as well the IP address of the platform, an example of which is shown in Figure 5c.

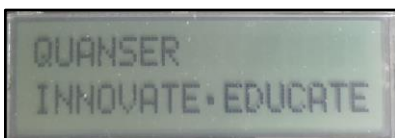

a. First bootup message



b. Second bootup message




c. Default display message

Figure 5. LCD startup and default messages


## E. Turning OFF the QCar

To turn OFF the QCar, firmly but gently press on the red power button once. This will display the message shown in Figure 6a. **This is the recommended process to initiate shutdown normally.** Please wait for the QCar to fully shut down before disconnecting the battery.

If for some reason the QCar will not shut down by clicking the power button, press and hold the power button for approximately 1 second. The message in Figure 6b will then be displayed. Upon releasing the power button, the power to the QCar computer will be cut off immediately and the LCD will go blank, equivalent to a hard shutdown. This is **not recommended** and should not be required during normal operation.

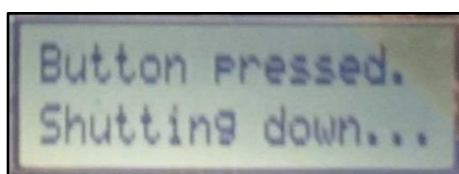

a. Normal shutdown message

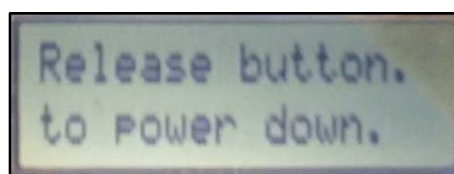

b. Forceful shutdown message

Figure 6. LCD shutdown message

## F. Low-battery and auto-shutdown

Lithium Polymer (LiPo) batteries can be damaged if discharged below a threshold voltage. A low battery warning will be displayed if the battery voltage is below **10.5 V**, as shown in Figure 7a. You should save your work and begin recharging the battery.

To prevent damage to the battery, the QCar will shut down automatically if the battery voltage gets below **10.0 V**, displaying the message shown in Figure 7b. It will first attempt a normal shutdown, but if it is unable to do so, it will disconnect the power after a short period.



a. Low battery warning message



b. Automatic shutdown

Figure 7. Low voltage protection and automatic shutdown

© Quanser Inc., All rights reserved.



Solutions for teaching and research. Made in Canada.